




## Digital High Temperature Modular Block

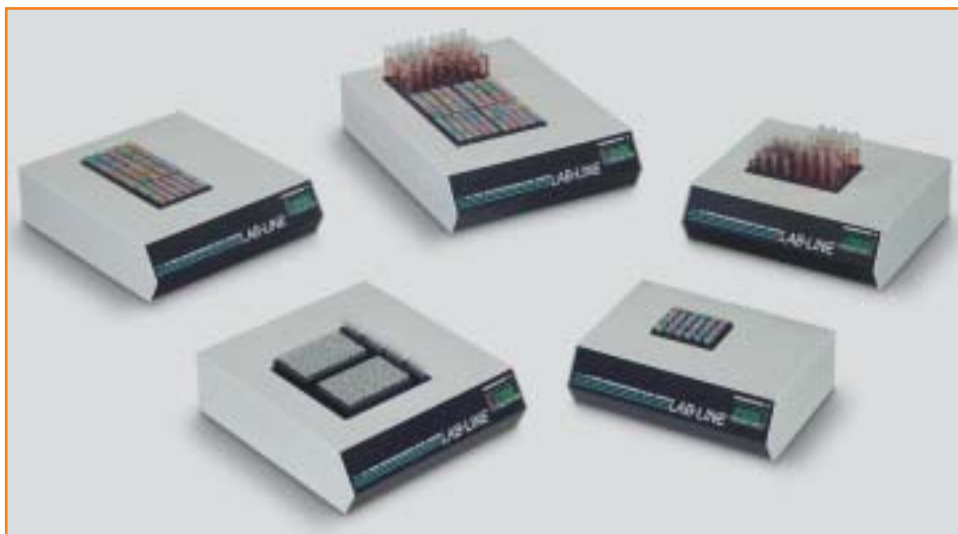

High Temperature Modular Blocks

**Includes FREE Built-in Temperature Probe**

### Temperature Range to 200°C!

- Choose from 1, 2, 3, 4, or 6 block configurations
- Temperature sensing probe system ensures accurate block temperature
- ☆ 24 Month Warranty

#### Product Description

- Provide uniform dry heating for cuvettes, microcentrifuge tubes, titer plates and test tubes.
- Close contact of tubes to block walls provides excellent heat transfer.
- Blocks include thermometer well for monitoring block temperature. Modular blocks sold separately.
- Solid black anodized aluminum blocks provide greater heat retention.
- Compact design conserves valuable bench top space.
- Powder coated steel body construction.
- Block puller included with all models.
- Use only Lab-Line modular blocks.

#### Operation

- Temperature uniformity  $\pm 0.2^\circ\text{C}$ .
- Temperature control of  $\pm 0.1^\circ\text{C}$ .

#### Applications

- Food sciences
- Hematology
- Incubation
- Chemistry
- Denaturing
- Pathology

### PRODUCT SPECIFICATIONS

Model #	Overall Dimensions Inches (cm)			Shipping Weight Lb. (kg)	Model #	Overall Dimensions Inches (cm)			Shipping Weight Lb. (kg)
	W	H	D			W	H	D	
2005	12.2 (31)	3.9 (10)	11.6 (19)	5 (2)	2008	12.2 (31)	3.9 (10)	15.9 (40)	8 (4)
2005-1	12.2 (31)	3.9 (10)	11.6 (19)	5 (2)	2008-1	12.2 (31)	3.9 (10)	15.9 (40)	8 (4)
2006	12.2 (31)	3.9 (10)	11.6 (19)	6 (3)	2009	12.2 (31)	3.9 (10)	19.6 (50)	10 (5)
2006-1	12.2 (31)	3.9 (10)	11.6 (19)	6 (3)	2009-1	12.2 (31)	3.9 (10)	19.6 (50)	10 (5)
2007	12.2 (31)	3.9 (10)	17.2 (44)	8 (4)	<i>All models include heavy duty 3-wire line cord and block puller; 120 volt models include plug.</i>				
2007-1	12.2 (31)	3.9 (10)	17.2 (44)	8 (4)					

### ORDERING INFORMATION

Model #	Volts	Hz	Amps	Watts	# of Modular Blocks Units Can Hold	US List Price	Model #	Volts	Hz	Amps	Watts	# of Modular Blocks Units Can Hold	US List Price
2005-1	240	50/60	0.8	200	1	\$933.00	2008-1	240	50/60	1.7	400	4	\$1,265.00
2006	120	50/60	2.1	250	2	\$1,037.00	2009	120	50/60	4.2	500	6	\$1,460.00
2006-1	240	50/60	1.0	250	2	\$1,037.00	2009-1	240	50/60	2.1	500	6	\$1,460.00
2007	120	50/60	2.5	300	3	\$1,163.00	<i>All models include heavy duty 3-wire line cord and block puller; 120 volt models include plug.</i>						
2007-1	240	50/60	1.3	300	3	\$1,163.00							