


Electrothermal Heating Tapes/Mats



Quartz to 800°C (1472°F)

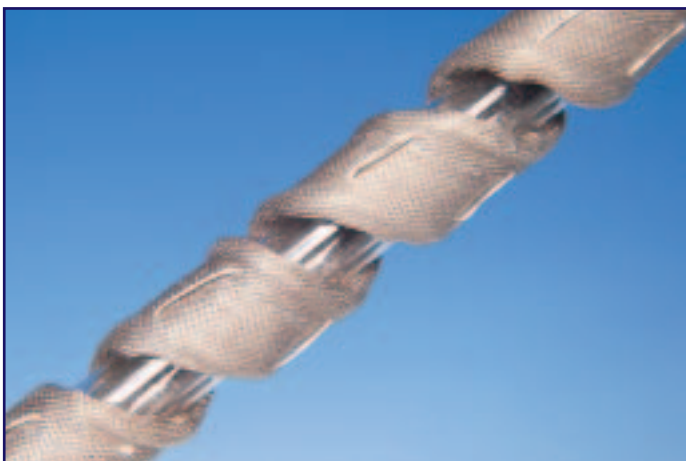

Model # HT75502*



*NOTE: All tapes are provided without controller. You must use a suitable controller or damage to your tape will occur. See controllers following heating mantle section. Tapes are not water resistant and should always be used in dry conditions and never for outdoor applications.
* To order 115V units add a suffix of X1 to the model number. i.e. HT75503 is 230V and HT75503X1 is 115V.*

Our Highest Temperature!

- Quartz Heat Tape
- Designed for many applications
- Temperature ranges to 800°C (1472°F)

HT7, HT9 Series - Up to 450°C (842°F)

- Constructed of high temperature resistant glass fabric for temperatures up to 450°C
- Inorganic core double glass covered heating element protected by 2 layers of glass fabric
- Grounded wire makes HT9 tapes safe for applying to metal pipes
- Free moving and self-adjusting elements
- Eyelets and tie-cord provide secure fixing

HT5 Series - Up to 800°C (1472°F)

- Quartz braided element is woven into an outer quartz fabric cover
- Recommended for high temperature (up to 800°C) applications
- Laces are provided at each end for secure fixing
- A 2 ft. (0.61m) lead with terminal block provided

HT6 Series - Up to 450°C (842°F)

- Glass fiber insulation
- A .5 meter cold lead with connector provide
- Safe for applying to metal pipes
- Laces are provided at each end for secure fixing
- Fixing to single and double curvature surfaces

ORDERING INFORMATION

Model #	Length ft. (m)	Watts	Volts	US List Price	Model #	Length ft. (m)	Watts	Volts	US List Price
HT7 & HT9 Series, Standard width 1" (25 mm) - Temperatures to 450°C (842°F)					HT5 Series, Standard width 1" (25 mm) - Temperatures to 800°C (1472°F)				
HT75502 ¹	2 (0.61)	100	115/230	\$87.00	HT550 ¹	2 (0.61)	200	115/230	\$276.00
HT75503	3 (0.91)	150	230	\$108.00	HT551 ¹	4 (1.22)	400	115/230	\$366.00
HT75504 ¹	4 (1.22)	200	230	\$117.00	HT552 ¹	6 (1.83)	600	115/230	\$487.00
HT75506 ¹	6 (1.83)	300	230	\$185.00	HT553 ¹	8 (2.44)	800	115/230	\$665.00
HT75508 ¹	8 (2.44)	400	230	\$193.00	HT558 ¹	16 (4.88)	1600	115/230	\$1,141.00
HT75512 ¹	12 (3.66)	600	230	\$228.00	HT6 Series, Standard width .5 (12.5) and 1" (25 mm) - Temperatures to 450°C (842°F)				
HT75516 ¹	16 (4.88)	800	230	\$286.00	HT640X1 ¹	.5 (12.5) x 2 (0.61)	100	115	\$129.00
HT75524 ¹	24 (7.32)	1200	230	\$393.00	HT641X1 ¹	.5 (12.5) x 4 (1.22)	200	115	\$143.00
HT75532 ¹	32 (9.76)	1600	230	\$667.00	HT642X1 ¹	.5 (12.5) x 6 (1.83)	300	115	\$154.00
HT95502 ¹	2 (0.61)	100	230	\$166.00	HT643X1 ¹	.5 (12.5) x 8 (2.44)	400	115	\$179.00
HT95503 ¹	3 (0.91)	150	230	\$179.00	HT650X1 ¹	1 (25) x 2 (0.61)	200	115	\$143.00
HT95504 ¹	4 (1.22)	200	230	\$215.00	HT651X1 ¹	1 (25) x 4 (1.22)	400	115	\$166.00
HT95506 ¹	6 (1.83)	300	230	\$228.00	HT652X1 ¹	1 (25) x 6 (1.83)	600	115	\$228.00
HT95508 ¹	8 (2.44)	400	230	\$251.00	HT653X1 ¹	1 (25) x 8 (2.44)	800	115	\$251.00
HT95512 ¹	12 (3.66)	600	230	\$337.00					
HT95516	16 (4.88)	800	230	\$397.00					
HT95524 ¹	24 (7.32)	1200	230	\$543.00					
HT95532 ¹	32 (9.76)	1600	230	\$665.00					

¹ CE Marked