


Thermolyne Heating Tapes/Mats



Silicone Rubber

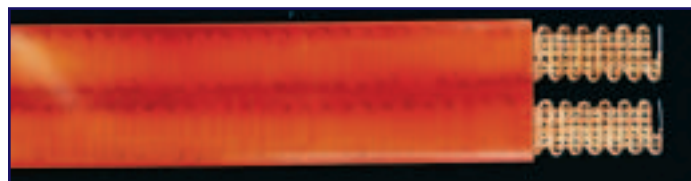

CE

Silicone Rubber Extruded Tapes

This tape is constructed of standard, insulated Fiberglass heating tapes completely encapsulated in void-free high strength, high-temperature silicon rubber. Designed to be moisture, chemical and abrasion resistant.

Grounded Silicone Rubber Extruded Tapes

These heavy insulated grounded heating tapes consist of a resistance wire wrapped in a layer of kapton tape, braided by two more layers of fiberglass insulation, then braided by a ground wire and covered with Silicone Rubber. Designed for strict grounding requirements, these tapes are moisture, chemical and abrasion resistant.

Built-in Controller, Grounded Tapes

These heating tapes have a built-in controller and are constructed of standard, insulated Fiberglass. Completely encapsulated in void-free high-strength, high-temperature silicon rubber, these tapes are designed to be moisture, chemical and abrasion resistant.

ORDERING INFORMATION

Model #	Dimensions Inches (cm x m)	Watts	Volts	US List Price
Silicone rubber 260°C (500°F) Standard Insulated				
BS0051-020	.5" x 2' (1.27 x 0.6)	52	120	\$71.00
BS0051-040	.5" x 4' (1.27 x 1.2)	104	120	\$78.00
BS0051-060	.5" x 6' (1.27 x 1.83)	156	120	\$92.00
BS0051-080	.5" x 8' (1.27 x 2.44)	209	120	\$102.00
BS0051-100	.5" x 10' (1.27 x 3.05)	261	120	\$116.00
BS0051-120	.5" x 12' (1.27 x 3.66)	313	120	\$131.00
BS0101-020	1" x 2' (2.5 x 0.6)	104	120	\$96.00
BS0101-040	1" x 4' (2.5 x 1.22)	209	120	\$118.00
BS0101-060	1" x 6' (2.5 x 1.83)	313	120	\$139.00
BS0101-080	1" x 8' (2.5 x 2.44)	418	120	\$162.00
BS0101-100	1" x 10' (2.5 x 3.05)	522	120	\$183.00
BS0101-120	1" x 12' (2.5 x 3.66)	627	120	\$211.00
BS0101-140	1" x 14' (2.5 x 4.3)	731	120	\$237.00
BS0101-160	1" x 16' (2.5 x 4.87)	836	120	\$258.00
BS0101-200	1" x 20' (2.5 x 6.1)	1045	120	\$296.00
BS0201-020	2" x 2' (5.1 x 0.6)	209	120	\$125.00
BS0201-040	2" x 4' (5.1 x 1.2)	418	120	\$177.00
BS0201-060	2" x 6' (5.1 x 1.83)	627	120	\$251.00
BS0201-080	2" x 8' (5.1 x 2.44)	836	120	\$309.00
BS0201-100	2" x 10' (5.1 x 3.05)	1045	120	\$342.00
BS0201-140	2" x 14' (5.1 x 4.3)	1463	120	\$391.00
BS0201-160	2" x 16' (5.1 x 4.87)	1672	120	\$446.00
BS0201-200	2" x 20' (5.1 x 6.1)	2090	120	\$546.00

Tape Tie Cords not included.

Model #	Volts	Watts	Dimensions Inches (cm x m)	Range	US List Price
Grounded Extruded Standard Insulated Tapes - Temperatures to 232°C (450°F)					
BS0101-020LG	120	104	1" x 2' (2.5 x 0.6)	232°C	\$231.00
BS0101-040LG	120	209	1" x 4' (2.5 x 1.2)	232°C	\$303.00
BS0101-060LG	120	313	1" x 6' (2.5 x 1.83)	232°C	\$366.00
BS0101-080LG	120	418	1" x 8' (2.5 x 2.44)	232°C	\$433.00
BS0101-100LG	120	522	1" x 10' (2.5 x 3.05)	232°C	\$455.00
BS0102-020LG	240	104	1" x 2' (2.5 x 0.6)	232°C	\$231.00
BS0102-040LG	240	209	1" x 4' (2.5 x 1.22)	232°C	\$303.00
BS0102-060LG	240	313	1" x 6' (2.5 x 1.83)	232°C	\$367.00
BS0102-080LG	240	418	1" x 8' (2.5 x 2.44)	232°C	\$433.00
BS0102-100LG	240	522	1" x 10' (2.5 x 3.05)	232°C	\$455.00

HT 88 Series, Grounded Extruded Tapes - Temperatures to 200°C (392°F)

HT8804X1*	115	80	1" x 4' (2.5 x 1.22)	200°C	\$151.00
HT8806X1*	115	120	1" x 6' (2.5 x 1.83)	200°C	\$168.00
HT8808X1*	115	160	1" x 8' (2.5 x 2.44)	200°C	\$204.00
HT8812X1*	115	240	1" x 12' (2.5 x 3.66)	200°C	\$203.00
HT8816X1*	115	320	1" x 16' (2.5 x 4.87)	200°C	\$264.00
HT8824*	230	480	1" x 24' (2.5 x 7.31)	200°C	\$331.00

Built-in Controller, Grounded Heat Tapes - Temperatures to 260°C (500°F)

BSAT101-020	120		1" x 2' (2.5 x 0.6)	316°C	\$126.00
BSAT101-040	120		1" x 4' (2.5 x 1.2)	316°C	\$139.00
BSAT101-060	120		1" x 6' (2.5 x 1.83)	316°C	\$167.00
BSAT101-080	120		1" x 8' (2.5 x 2.44)	316°C	\$179.00
BSAT101-100	120		1" x 10' (2.5 x 3.05)	316°C	\$190.00

Grounded LG Tapes supplied without plugs.

* Tapes are also available in 230 V and are CE marked.