


## Lab-Line Ultrasonic Cleaners



# Digital Aquawave™

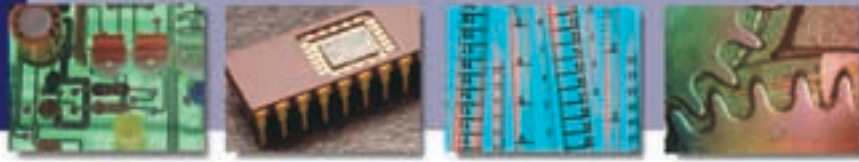

Model 9331

## Digital Models With More Options!

- LED digital displays for temperature and time
- 40 KHz cleaning power
- Stainless steel construction
- Available in sizes from 2.75 l to 28 l
- ☆ 24 Month Warranty

### Product Description

- Combine efficiency and thoroughness of ultrasonic cleaning with the accuracy and ease-of-use of digital technology.
- Durable stainless steel tank and cabinet.
- 40 kHz operation for optimal cleaning results.
- Sandwich-type transducers tolerate all-day continuous use.
- Adjustable ultrasound power level—up to 140% due to second transducer.
- Digital temperature and timer display.
- Non-volatile memory maintains all settings after shut-down and in the event of power loss.
- Adjustable heat—from 20°C to 80°C in 1.0°C increments.

- Degas function—increases medium cleansing action and shortens total duration of cleaning cycle.
- Quick-set control knobs.
- Overheat protection thermostat.
- Built-in timers can be set for timed cycles of up to 90 minutes or continuous operation.
- Drains included on large models.

### Applications

- Lab glassware
- Surgical and dental instruments
- Motors
- Scientific instruments
- PC boards

### PRODUCT SPECIFICATIONS

Model #	Tank Dimensions Inches (cm)			Overall Dimensions Inches (cm)			Model #	Tank Dimensions Inches (cm)			Overall Dimensions Inches (cm)		
	W	H	D	W	H	D		W	H	D	W	H	D
9331	9.4 (23.9)	3.9 (10.0)	5.4 (13.7)	10.4 (26.4)	8.5 (21.6)	6.9 (17.5)	9333	11.8 (30.0)	7.9 (20.1)	9.4 (23.9)	17.7 (45.0)	15.7 (39.9)	11.4 (29.0)
9331-1	9.4 (23.9)	3.9 (10.0)	5.4 (13.7)	10.4 (26.4)	8.5 (21.6)	6.9 (17.5)	9333-1	11.8 (30.0)	7.9 (20.1)	9.4 (23.9)	17.7 (45.0)	15.7 (39.9)	11.4 (29.0)
9322	19.9 (50.5)	3.9 (10.0)	5.4 (13.7)	22.2 (56.5)	8.5 (21.5)	7.9 (20.0)	9334	12.9 (32.7)	7.9 (20.1)	11.8 (30.0)	17.7 (45.0)	15.7 (39.9)	15.4 (39.4)
9322-1	19.9 (50.5)	3.9 (10.0)	5.4 (13.7)	22.2 (56.5)	8.5 (21.5)	7.9 (20.0)	9334-1	12.9 (32.7)	7.9 (20.1)	11.8 (30.0)	17.7 (45.0)	15.7 (39.9)	15.4 (39.4)
9323	11.8 (30.0)	5.9 (15.0)	9.4 (24.0)	15.5 (39.4)	14.7 (37.2)	10.8 (27.5)	9335	19.9 (50.5)	7.9 (20.1)	11.8 (30.0)	23.9 (60.7)	15.7 (39.9)	15.9 (40.4)
9323-1	11.8 (30.0)	5.9 (15.0)	9.4 (24.0)	15.5 (39.4)	14.7 (37.2)	10.8 (27.5)	9335-1	19.9 (50.5)	7.9 (20.1)	11.8 (30.0)	23.9 (60.7)	15.7 (39.9)	15.9 (40.4)

### ORDERING INFORMATION

Model #	Capacity Liters/Gal	Drain	Electrical (50/60 Hz)		US List Price	Model #	Capacity Liters/Gal	Drain	Electrical (50/60 Hz)		US List Price
			Volts	Watts					Volts	Watts	
9331	2.8/7.4	No	120	300	\$626.00	9333	12.9/3.4	Yes	120	600	\$1,464.00
9331-1	2.8/7.4	No	240	300	\$626.00	9333-1	12.9/3.4	Yes	240	600	\$1,464.00
9322	6.8/1.8	Yes	120	360	\$1,125.00	9334	18.0/4.8	Yes	120	600	\$1,810.00
9322-1	6.8/1.8	Yes	240	360	\$1,137.00	9334-1	18.0/4.8	Yes	240	600	\$1,810.00
9323	9.5/2.5	Yes	120	360	\$1,160.00	9335	28.0/7.4	Yes	120	600	\$2,000.00
9323-1	9.5/2.5	Yes	240	360	\$1,160.00	9335-1	28.0/7.4	Yes	240	600	\$2,000.00