



Flow Indicator

Provides Immediate Visual Indication of Water, Air or Gas Flow

- Highly visible red nylon ball
- Ball rotates in a circular channel during flow
- Continuity and relative flow rate indicated by speed of ball's rotation

The minimum flow rate to actuate the ball with water is 50ml per minute in a horizontal position and 140ml per minute in a vertical position. With air, the minimum flow rate is 3 liters per minute in a horizontal position and 8 liters per minute in a vertical position. Maximum working pressure is 10 p.s.i.g. (69 kN/m²). The polycarbonate indicator measures 39 x 57mm (1½ x 2¼") with two 23mm (.09") long tube extensions to accept 6.4mm (¼") to 9.5mm (¾") I.D. tubing. An arrow shows direction of flow. **Per each, 12 per case.**



CATALOG NO.	PRICE/EACH
H19935-0000	\$14.85



Flow Monitor

Red Paddle Wheel Clearly Indicates Liquid or Gas Flow

- Clear polystyrene outer case reveals bright red paddle wheel
- Paddle wheel spins freely with flow

This superior in-line device indicates liquid or gas flow from 95ml/sec. down to 5ml/sec. (or down to

3ml/sec. in fluid with air bubbles removed). Not for use with corrosive fluids or over 64°C (150°F). Dimensions are 90 x 16 x 55mm (3½ x ¾ x 2¼"). Serrated tubulations accept 10 to 14mm (¾ to ½") ID tubing. **Per each, 12 per case.**



CATALOG NO.	PRICE/EACH
H19937-0000	\$35.25



Roto-Flo™ Flow Indicators

Colored Paddle wheels Go with the Flow

- Provide immediate visualization of flow in a tubing line
- Brightly colored polypropylene paddle wheel spins when gas or liquid flows through a system


Body molded of clear polystyrene with smooth connecting barbs for easy installation. Internal flow channel diameter is 5mm, paddle wheel shell is 25mm, overall length is 51mm. Maximum recommended pressure is 10 p.s.i.g. Minimum flow: 0.5 LPM. Maximum flow: 10 LPM. Not for use with corrosive chemicals or fluids with a temperature in excess of 64°C (150°F). **Per each, 12 per case.**



CATALOG NO.	FITS TUBING I.D.	TUBE ENDS	WHEEL COLOR	PRICE/EACH
H19937-0001	4mm to 7mm (¼ to ½")	straight	Blue	\$13.20
H19937-0002	4mm to 7mm (¼ to ½")	Barbed	Red	13.20
H19937-0003	7mm to 11mm (¼ to ½")	Barbed	Green	13.20

1-800-4BEL-ART

FUME HOODS

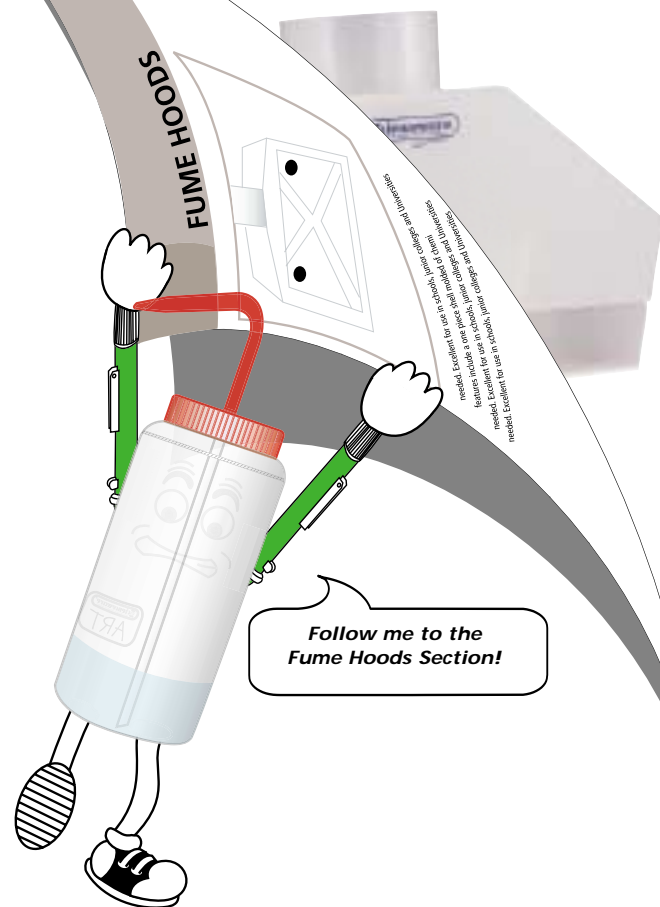

**TUBING,
FITTINGS
and
FAUCETS,
See Pages,
163-165** →