

### Pipette Jars (A)

For Cleaning and Storage

Leakproof, molded polyethylene Pipette Jars store and soak pipettes and will withstand all common cleaning solutions at room temperature, including sulfuric acid/dichromate. The Scienceware® H17106-0000 or H17107-0000 Pipette Basket (below), filled with pipettes, fits into the corresponding jar. 152mm (6") in diameter. **Per each, 4 per case.**



CATALOG NO.	HEIGHT	PRICE/EACH
F17130-0018	46cm (18")	\$61.80
F17130-0024	61cm (24")	64.90

### Pipette Baskets (B)

Avoid Costly Breakage

Cylindrical pipette baskets have no handle to obstruct the insertion and retrieval of pipettes. Instead, the sides have cut-out hand grips for easy handling. Gusher holes provide quick drainage when lifted from a soaking jar, and a foam disk on the bottom cushions pipette tips to prevent breakage. Polyethylene baskets are 102mm (4") I.D. **Per each, 6 per case.**



CATALOG NO.	FOR PIPETTES	HEIGHT	PRICE/EACH
H17106-0000	46cm (18")	58cm (23")	\$75.00
H17107-0000	61cm (24")	81cm (32")	81.30

### Pipette Rinsers (C)

Effective, Reliable Design Ensures Safe, Efficient Operation

High level water entry provides an air gap between the maximum water level and the end of the water inlet tube ensuring safe and efficient operation. The Scienceware® polyethylene pipette rinser provides a simple, automatic method of rinsing large numbers of pipettes at one time. The "tube within a tube" siphoning system assures rapid, complete and automatic siphoning and refilling over a wide range of inlet water pressures and flow rates. The design minimizes the danger of siphon stalling at certain water inflow rates. It prevents any possible backflow into the potable water supply in the event of a pressure interruption. Each rinser includes a 305mm (12') long piece of 2.5cm (1") I.D. drain tube and 1.5M (5ft) of 12.5mm (½") I.D. flexible water inlet tube. Octagonal base is 26cm (10½") across. Cylinder is 152mm (6") I.D., slant top, accepts Scienceware® H17106 & H17107 pipette baskets.



CATALOG NO.	FOR PIPETTES	HEIGHT	PRICE/EACH
F17131-0018	46cm (18")	64cm (25½")	\$175.00
F17131-0024	61cm (24")	79cm (31½")	179.00




### Cleanware™ Stainless Steel Sponge

Cleans without Scratching

Scrub glassware and stainless steel products sparkling clean without damage to labware. For use with or without detergents. There are no sharp edges to scratch surfaces or irritate hands. Measuring roughly 64 x 100mm (2.5 x 4"), these pliant sponges are flexible enough to clean hard to reach areas thoroughly. Grade 430 steel wool. **2 per bag, 24 bags per case.**



CATALOG NO.	PRICE/BAG
F17085-0000	\$12.00

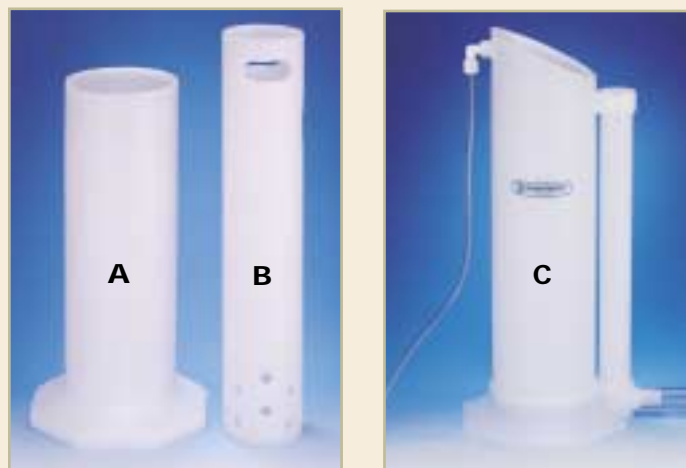

### The Complete Cleanware™ Pipette Rinsing Systems

The TRIO

The 46cm (18") TRIO for pipettes up to 46cm (18") long. Includes one each of items A (F17130-0018), B (H17106-0000), and C (F17131-0018) shown above.



CATALOG NO.	PRICE/EACH
H17121-0000	\$264.99

THE 61cm (24") TRIO for pipettes up to 61cm (24") long. Includes one each of items A (F17130-0024), B (H17107-0000), and C (F17131-0024) shown above.



CATALOG NO.	PRICE/EACH
H17122-0000	\$295.20

Choose from a Variety of Pipette Storage Racks and Trays, See Page 97-98

### Test Tube Baskets

For Rinsing, Transport and Storage

Ideal for rinsing glassware, these polypropylene baskets with lids are simply inverted for easy draining. They can also be used for transporting small laboratory animals. Steam autoclavable at 121°C (250°F). **Per each, 6 per case.**



CATALOG NO.	DIMENSIONS (LENGTH X WIDTH X HEIGHT)	PRICE/EACH
F18737-0000	127 x 102 x 102mm (5 x 4 x 4")	\$11.85
F18737-0010	With lid. 127 x 102 x 102mm (5 x 4 x 4")	14.45
F18738-0000	152 x 152 x 152mm (6 x 6 x 6")	15.45
F18738-0010	With lid. 152 x 152 x 152mm (6 x 6 x 6")	16.75
F18739-0000	229 x 229 x 229mm (9 x 9 x 9")	18.25
F18739-0010	With lid. 229 x 229 x 229mm (9 x 9 x 9")	20.60

### The Glascribe® Pen

For Marking Laboratory Glassware

The Glascribe® Pen, with retractable tungsten carbide tip, lets you inscribe on glass, ceramics, and plastic. The pen has a pocket clip so it is always handy when needed. **Per each, 12 per case.**



CATALOG NO.	PRICE/EACH
H44150-0000	\$17.25

See the Labeling Section for Felt Tip Glass Marking Pens, Pages 69-71