


Pipette Pump™ III Pipettors

Patent Pending

The ergonomic design of every Pipette Pump™ III Pipettor enhances comfort, utility and dispensing precision while helping to reduce the risk of injury from repetitive and awkward pipetting motions. The angle of the inserted pipette and its placement in front of the hand allows for use with the arm in a more natural, lowered position unlike other pipettors, which require arm and elbow elevation for many pipetting actions. The compact design adds only 152 mm (6") to the pipette length allowing use in fume hoods and controlled environment glove boxes. The pipettors are lightweight for one-handed operation and minimal strain on hand and wrist.

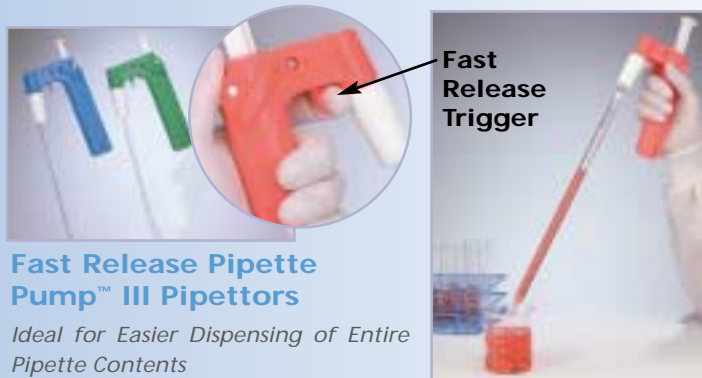

Fast Release Pipette Pump™ III Pipettors

Ideal for Easier Dispensing of Entire Pipette Contents

The lightweight Fast Release Pipette Pump™ III Pipettor is designed for fast and accurate one-handed operation. A fast release trigger dispenses the entire pipette contents. An optimally located thumb wheel rotates easily for precision aspirating or dispensing, down to a volume of one drop at a time – all while using a more relaxed grip. The thumb wheel also can be used to “blow out” the pipette. Available in three different pipette volumes, 2, 10 and 25 ml, each color coded for easy identification. A self-stabilizing collar inside the silicone chuck holds either glass or plastic pipettes. Pipettors are resistant to acids, alkalis and most solvents and can be disassembled for cleaning. Overall length 171 mm (6.75"). **Per each, 12 per case.**

CATALOG NO.	FOR PIPETTES UP TO	COLOR	PRICE/EACH
F37904-0002	2ml	blue	\$17.25
F37904-0010	10ml	green	17.50
F37904-0025	25ml	red	23.50

Economy Pipette Pump™ III Pipettors

Ergonomic Design – Affordable Price !

The ergonomic design of the Economy Pipette Pump™ III Pipettor enhances comfort, utility and dispensing precision while helping to reduce the risk of injury from repetitive and awkward pipetting motions. This lightweight, affordable device provides accurate, trouble-free pipetting with simple, one-handed operation. The optimally located thumb wheel rotates easily for precision aspirating or dispensing, down to a volume of one drop at a time – all while using a more relaxed grip. The thumb wheel also can be used to “blow out” the pipette. Available in three different pipette volumes, 2, 10 and 25 ml, each color coded for easy identification. It resists acids, alkalis and most solvents and can be easily disassembled for cleaning. Overall length 171 mm (6.75"). **Per each, 12 per case.**

CATALOG NO.	FOR PIPETTES UP TO	COLOR	PRICE/EACH
F37902-0002	2ml	blue	\$15.75
F37902-0010	10ml	green	16.50
F37902-0025	25ml	red	22.50

NEW

Fast Release Pipette Pump™ II Pipettors

Features Rapid Release Valve for One-Handed Operation

Designed for fast and efficient single -handed operation with glass or plastic pipettes. Pipettes are easily inserted into a unique collar with stabilizing fingers to hold pipette securely. Rotating the knurled thumb wheel draws liquid to the desired volume. To empty rapidly, apply light pressure to the fast release lever. For gradual dispensing, rotate the thumb wheel. Easily disassembled for cleaning. Sizes are color coded. **Per each, 12 per case.**

CATALOG NO.	FOR PIPETTES UP TO	COLOR	PRICE/EACH
F37911-1002	2ml	blue	\$17.25
F37911-1010	10ml	green	17.50
F37911-1025	25ml	red	23.50

Up to 25ml

Up to 10ml

Up to 2ml

Up to 0.2ml

Pipette Pump™ Pipettors

Color Coded by Size

This handy device provides safe, accurate and trouble free pipetting with simple, one-handed operation. Pipettes fit smoothly into its flexible tapered chuck. The thumbwheel on the side is rotated for precision filling or dispensing, and the plunger may be pressed for quick emptying. It resists acids and alkalis, and is easily disassembled for cleaning. Color of barrel denotes pipette size. **Per each, 12 per case.**

CATALOG NO.	FOR PIPETTES UP TO	COLOR	PRICE/EACH
F37896-0000	0.2ml	yellow	\$17.25
F37897-0000	2ml	blue	15.75
F37898-0000	10ml	green	16.50
F37899-0000	25ml	red	22.50

Spare parts for the Pipette Pump™ Pipettors, silicone rubber chuck and washers

FOR MODEL	PIPETTE CHUCK	PRICE/EACH	WASHERS (SET)	PRICE/EACH
F37896-0000	H93789-0060	\$1.80	H93789-0061 (sleeve)	\$1.59
F37897-0000	H93789-0070	3.00	H93789-0071 3/Set	1.35
F37898-0000	H93789-0080	3.09	H93789-0081 3/Set	1.95
F37899-0000	H93789-0090	3.30	H93789-0091 3/Set	2.19

Volac™ Universal Pipette Controllers

U.S. Patent No. 3,809,297

A Quick Repetitive Pipettor for Pre-Set Volumes

A superior device for quick, one-handed, repetitive pipetting of pre-set volumes. It accepts standard pipettes from 0.1 to 20ml with diameters of 4.5 to 8.7mm. The Volac™ Pipettor incorporates plunger action and a fast release pump with a fully enclosed mechanism for adjusting the volume setting. The unique free flow control permits fractional delivery. Color coded for safety and convenience, the entire unit is autoclavable. Spare parts included. **Per each, 6 per case.**

CATALOG NO.	DESCRIPTION	COLOR	PRICE/EACH
H37892-0000	For 0.1 to 2.0ml pipettes	red	\$154.50
H37893-0000	For 2.0 to 20 ml pipettes	blue	154.50

