

Adapt-A-Port® T Stopper Kits, Teflon® PTFE


The Safe-Lab® Adapt-A-Port® Kit is a unique adapter set consisting of three Teflon® PTFE pieces. The Teflon® PTFE body is threaded on the inside top to accept either the threaded vacuum port adapter or thermometer port adapter. Total interchangeability is achieved by simply unscrewing one threaded port adapter and replacing it with the other. The kit comes complete with O-rings. Not Autoclavable.



The kit contains:

(A) One standard taper Teflon® PTFE Adapt-A-Port® body with a Viton® O-ring for air-tight sealing. (B) One Teflon® PTFE threaded vacuum port adapter with a side angled serrated Teflon® PTFE nozzle for the tubing connections. (C) One Teflon® PTFE threaded thermometer port adapter with locking compression cap on top for straight stem thermometers or glass tubing. An O-ring on the bottom threaded section guarantees air-tight sealing. Items B and C are available only in the Adapt-A-Port® Kit. Not Autoclavable.

CATALOG NO.	JOINT SIZE	PRICE/EACH
F20304-1560	24/40	\$79.00
F20304-1565	29/42	82.00

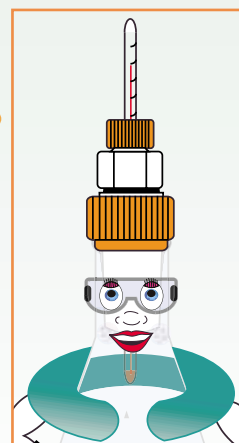

Safe-Lab® Adapt-A-Port® T Body

Teflon® PTFE Adapt-A-Port® body threaded on top to accept the Scienceware® F20306 threaded port adapters. Comes with polypropylene extracting nut. Not autoclavable.



CATALOG NO.	JOINT SIZE	PRICE/EACH
F20303-1635	24/40	\$37.50
F20303-1640	29/42	46.00



Safe-Lab® Thermometer T Joint Adapters

Teflon® PTFE body with Viton® O-ring for air tight vacuum work. The adapter has a threaded polypropylene screw cap on top with an inside O-ring for depth adjustment of 6-7mm for straight stem thermometers, glass tubing and other items. Rotating the extracting nut removes the adapter from flask. The size 10/18 adapter has no body O-ring. Not autoclavable.



CATALOG NO.	JOINT SIZE	PRICE/EACH
F20309-1675	10/18	\$30.00
F20309-1680	14/20	33.75
F20309-1710	19/22	38.50
F20309-1716	24/40	42.75



Safe-Lab® Therm-O-Vac® T Joint Adapter

Teflon® PTFE adapter with Viton® O-ring on the body, screw cap and hose connection. Not autoclavable.



CATALOG NO.	JOINT SIZE	HOLE OPENING	PRICE/EACH
F20307-1825	24/40	6mm	\$65.00



Safe-Lab® Joint Tubing Adapters

Teflon® PTFE body with Viton® O-ring. Polypropylene screw cap on top with inside O-ring for depth adjustment of glass tubing. Bottom inner joint has a threaded extracting nut for easy removal of adapter. Not autoclavable.



CATALOG NO.	JOINT SIZE	HOLE OPENING	PRICE/EACH
F20311-1717	24/40	8mm	\$47.00
F20311-1718	24/40	10mm	47.00



Safe-Lab® Therm-O-Vac® Port Adapters

Threaded Therm-O-Vac® Port Adapter with serrated side Teflon® PTFE nozzle. For use with Scienceware® F20303 Adapt-A-Port® Teflon® PTFE body. Not autoclavable.



CATALOG NO.	HOLE OPENING	PRICE/EACH
F20306-1610	6mm	\$42.75
F20306-1615	8mm	42.75