

# SELECTION GUIDE

## *Histology*



**OUR LARGE SIZE ECONOMICAL SLIDE-WARMER!**



**OUR COMPACT ECONOMICAL SLIDE-WARMER!**



**TOP OF THE LINE FLOAT BATH!**



**THE HOTTEST WAX DISPENSER IN THE MARKET!**

### SLIDE WARMERS

	Type	Operating Temp Range °C	W	Overall Dimensions Inches (cm) H	D	Page
<i>OUR LARGE SIZE ECONOMICAL SLIDE-WARMER!</i>	26000	Ambient to 65°C	14 (36) 25 (64) 25 (64)	14 (36) 8 (20) 8 (20)	3.5 (9) 3.5 (9) 4.25 (11)	122
<i>HIGHER TEMPERATURE SLIDE-WARMER WITH DRYING BARS!</i>	MH6600	Ambient to 100°C	15.6 (39.6)	7 (17.8)	4 (10.2)	123
<i>OUR COMPACT ECONOMICAL SLIDE-WARMER!</i>	26000	45° to 50°C	13 (33)	4 (10.1)	4 (10.1)	124

### FLOAT BATH

	Type	Operating Temp Range °C	Chamber Capacity Liters	Page
<i>LOW-PROFILE ECONOMICAL FLOAT BATH!</i>	26100	Ambient to 70°C	2.3	125
<i>TOP OF THE LINE FLOAT BATH!</i>	MH8500	Ambient to 75°C 40° to 80°C	2.5	126

### WAX DISPENSER

	Type	Operating Temp Range °C	Chamber Capacity Liters	Page
<i>THE HOTTEST WAX DISPENSER IN THE MARKET!</i>	MH8523	Ambient to 75°C	45	127

# Histology

## SLIDE WARMERS

Lab-Line

Model 26025



### OUR LARGE SIZE ECONOMICAL SLIDE-WARMER!

- — Temperature range from slightly above ambient to 65°C
- — Two sizes—accommodates up to 66 slides
- — Large slide capacity with space-saving design
- ★ — 24 month warranty



### PRODUCT DESCRIPTION

- Uniform heat transfer over entire heating surface.
- Anodized black heavy gauge aluminum work top resists corrosion and provides high-contrast background for viewing specimen.
- Powder coated heavy-duty steel body construction is durable and easy to clean.
- Thermostat-controller regulates temperature to 65°C.
- Model 26025 features LED display and cover.
- Model 26007 and Model 26022 are CE certified.
- Models 26005 and 26020 are ETLmarked.
- Models 26005 & 26007 have a slide capacity of sixty-five 1" x 3" (2.54 cm x 7.6 cm) or thirty-two 2" x 3" (5.1 cm X 7.6 cm).
- Models 26020, 26022, 26025 & 26025-1 have a slide capacity of sixty-six 1" x 3" (2.54 x 7.6 cm) or thirty-three 2" x 3" (5.1 x 7.6 cm).

### APPLICATION

- Designed for cytology, histology, and other biological applications in industrial and clinical lab environments.

### ORDERING INFORMATION

Model #	Electrical (50/60 Hz)		Overall Dimensions Inches (cm)			US List Price
	Volts	Watts	W	H	D	
26005	120	160	14 (36)	14 (36)	3.5 (9)	\$442.00
26007	240	160	14 (36)	14 (36)	3.5 (9)	\$442.00
26020	120	160	25 (64)	8 (20)	3.5 (9)	\$489.00
26022	240	160	25 (64)	8 (20)	3.5 (9)	\$489.00
26025	120	160	25 (64)	8 (20)	4.25 (11)	\$655.00
26025-1	240	160	25 (64)	8 (20)	4.25 (11)	\$655.00

### ACCESSORIES

Model #	Description	US List Price
26523	Metal dust cover for models 26022 and 26020	\$58.00

Model 26005



# Histology

## SLIDE DRYING BENCH

Electrothermal

MH6616



### HIGHER TEMPERATURE SLIDE-WARMER WITH DRYING BARS!

- — Temperature range from slightly above ambient to 100°C
- — Accepts up to 50 slides
- — Dry slides across the drying bars, angled on the bars or lying flat without bars

CE

### PRODUCT DESCRIPTION

- Painted aluminum top reduces glare, offering a high-contrast background for viewing specimen.
- Aluminum housing, black and white exterior is durable and easy to clean.
- Clear light indicates when power is on; amber light indicates when the heater is on.
- Temperature control is provided by a built-in energy regulation with a temperature range to 100°C.
- Reference dial and knob allow the temperature to be closely adjusted and regulated.
- Slide capacity: 50 (1" x 3") slides.
- Easily replaceable heating element.

### ORDERING INFORMATION


Model #	Electrical (50/60 Hertz)		Overall Dimensions Inches (cm)			Shipping Weight Lb. (kg)	US List Price
	Volts	Watts	W	H	D		
MH6616X1	115	150	15.5 (39)	3.5 (9)	7.12 (18)	5.3 (2.4)	\$743.00
MH6616	230	150	15.5 (3.9)	3.5 (9)	7.12 (18)	5.3 (2.4)	\$743.00

# Histology

## RH VIEW BOX

Thermolyne

Type LL16100

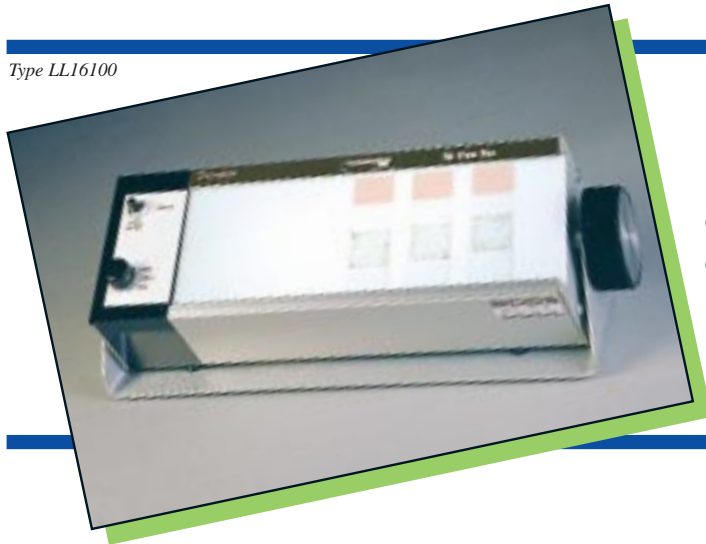

**OUR COMPACT ECONOMICAL SLIDE-WARMER!**

- — Soft, fluorescent, glare-free light
- — Built-in temperature indicator

### PRODUCT DESCRIPTION

- Built-in temperature indicator for easy, accurate monitoring of viewing area.
- Viewing area temperature adjusts easily to compensate for ambient temperature changes. (Range: 45°C to 50°C).
- Rh slide reaches temperature of 37° to 39°C in two minutes, eliminating wasted heat up time.
- Soft, fluorescent, glare-free light provides uniform illumination of slide contents.
- View box can be gently rocked using large knob on support cradle.
- Sturdy, easy-to-clean extruded aluminum case is compact; utilizes minimal bench space.
- Meets Rh slide test standards set by the AABB.
- Unit conveniently includes three-wire cord and plug.

### APPLICATIONS

- Blood typing.
- Rh determinations.
- Warming slides for gram-stain.
- Tissue typing.
- Reading enzyme antibody screens in microplate.

### PRODUCT SPECIFICATIONS

Type	Viewing Area Dimensions Inches (cm)		Overall Dimensions Inches (cm)			Shipping Weight Lb. (kg)
	W	D	W	H	D	
LL16100	9.50 (24.1)	3.06 (7.77)	13 (33)	4 (10.1)	4 (10.1)	4.25 (1.92)

### ORDERING INFORMATION

Model #	Volts	Electrical (50/60 Hertz) Amps	Watts	Viewing Area Operating Temp Range °C	US List Price
LL16115	120	0.24	30	45-50°C	\$349.00
LL16110-26	220	0.12	25	45-50°C	\$349.00

*-26 Models Supplied with European cord set.*

48 Hour or Sooner  
Express Shipping  
Guaranteed!\*



\*see inside front cover for details

# Histology

## TISSUE FLOAT BATH

Lab-Line

Model 26103



### LOW-PROFILE ECONOMICAL FLOAT BATH!

- — Controlled temperature up to 70°C
- — Seamless 2.3 liter aluminum chamber
- ★ — 24 month warranty

### PRODUCT DESCRIPTION

- Designed to thoroughly soften embedded tissue specimens without creating wrinkles, folds or distortions.
- One-piece aluminum chamber is easy to clean and won't leak into electrical component areas.
- 2.3 liter black finished interior reduces glare, improving visibility of tissue specimen.
- Top ledge of bath allows space for preheating or for evaporating away any remaining water on the glass slides.
- Chamber is completely encircled by sheathed heater for uniform, broad surface heating.
- Fiberglass-insulated body prevents heat loss and temperature fluctuations.
- Chamber capacity of 2.3 liters.
- Model 26106 has a LED temperature display and a hydraulic temperature control.
- Model 26106 has a removable 1.9 liter pyrex glass pan for easier cleaning.
- Model 26106 has a 6 watt fluorescent lamp that illuminates the bath pan to increase visibility.

### OPERATION

- Front mounted on/off switch and temperature control knob allow for complete control of the unit.
- Clear light indicates "Power On" and amber light indicates "Heater On."

### PRODUCT SPECIFICATIONS



Model #	Chamber Capacity (l)	Chamber Dimensions Inches(cm)		Overall Dimensions Inches(cm)		Shipping Weight Lb. (kg)
		H	Diameter	H	Diameter	
26103	2.3	2.5 (6)	8.5 (22)	3.5 (9)	13.25 (34)	9.5 (4.3)
26104	2.3	2.5 (6)	8.5 (22)	3.5 (9)	13.25 (34)	9.5 (4.3)

Model #	Internal Dimensions Inches (cm)			Overall Dimensions Inches (cm)			Shipping Weight Lb. (kg)
	W	H	D	W	H	D	
26106	10.75 (27.3)	7 (17.8)	1.75 (4.5)	20 (50.8)	23 (58.4)	15 (38.1)	18 (8.2)

### ORDERING INFORMATION

Model #	Electrical (50/60 Hz)		US List Price
	Volts	Watts	
26103	120	300	\$376.00
26104	240	300	\$376.00
26106	120	600	\$750.00

Model 26106



# Histology

## PARAFFIN SECTION MOUNTING BATH

**Electrothermal**

MH8515



### TOP OF THE LINE FLOAT BATH!

- — Comfortable and easy manipulation of sections
- — Constant clear visibility
- — Coated aluminum housing

CE

### PRODUCT DESCRIPTION

- Easy manipulation of the section avoids problems with wrinkles, folds and distortions.
- 2.5 liter one piece black interior is easy to clean and provides constant clear visibility of the sample.
- Aluminum exterior with polypropylene lid is durable and easy to clean.
- Easily replaceable heating element.
- Built-in controller regulates the temperature range of ambient to 75°C (Model MH 8513B); 40°C to 80°C (Model MH8514B).

### OPERATION

- Front mounted on/off switch and temperature control knob allows for complete control of the unit.
- Clear light indicates "Power On." Red light warns of over-temperature. Amber light indicates "Heater On."
- Excess or damaged surfaces can be easily removed by skimming a filter paper over the water surface.

### SAFETY

- Built-in safety device shuts off the energy supply to the heater in the event of an overheat condition such as boiling dry.

### PRODUCT SPECIFICATIONS

Model #	Chamber Capacity (l)	Chamber Dimensions Inches (cm)		Overall Dimensions Inches (cm)		Shipping Weight Lb. (kg)
		H	D	H	D	
MH8515X1	2.25	3.0 (7)	8.7 (22)	5 (13)	13 (33)	6 (2.7)
MH8515	2.25	3.0 (7)	8.7 (22)	5 (13)	13 (33)	6 (2.7)

### ORDERING INFORMATION

Model #	Electrical (50/60 Hertz)		US List Price
	Volts	Watts	
MH8515X1	115	240	\$929.25
MH8515	230	240	\$929.25

# Histology

## PARAFFIN WAX DISPENSER

Electrothermal

MH8523BX1



### THE HOTTEST WAX DISPENSER IN THE MARKET!

- — Provides immediate "Wax on Tap"
- — For use in Pathology laboratories, etc.
- — Thermally insulated tap heater with independent temperature control eliminates wax solidification



### PRODUCT DESCRIPTION

- Anodized black interior, stove enameled black and white aluminum exterior and a black metal lid are durable and easy to clean.
- Capacity of 1 gallon (45 liters).
- Built-in controller provides temperature control from ambient to 75°C.
- A non-drip swivel lever with replaceable filter makes the dispensing of wax easy.

### OPERATION

- The clear light indicates "Power On." The red light warns of over-temperature. Amber neon lights indicate "Bath Heater On" and "Tap Heater On."
- Tap heater maintains temperature within the tap to keep the wax molten.

### SAFETY

- Safety over-temperature at 105°C.
- Non-resettable thermal fuse.

### PRODUCT SPECIFICATIONS

Model #	Chamber Capacity (l)	Chamber Dimensions Inches (cm)		Overall Dimensions Inches (cm)		Shipping Weight Lb. (kg)
		H	Diameter	H	Diameter	
MH8523BX1	(4.5)	8.25 (15)	8.25 (15)	12 (30)	11 (28)	10 (4.5)
MH8523B	(4.5)	8.25 (15)	8.25 (15)	12 (30)	11 (28)	10 (4.5)

### ORDERING INFORMATION

Model #	Electrical (50/60 Hertz)		US List Price
	Volts	Watts	
MH8523BX1	115	175	\$2079.00
MH8523B	230	175	\$2079.00