

NEW

Poxygrid® Imhoff Cone Supports

These lightweight Imhoff Cone Supports hold three (or four) settling cones in a single row for clearer viewing. Cone tips are held up off the benchtop for easier viewing and to allow plug removal if required. The wire frame is epoxy coated for durability and easy cleaning. A wide base provides stability without crowding the benchtop.

CATALOG NO.	DESCRIPTION	DIMENSIONS	PRICE/EACH
F38993-0003	3 place	44.5 x 17.1 x 40.6cm H (17½ x 6¾ x 16"H)	\$45.95
F38993-0004	4 place	57.7 x 17.1 x 40.6cm H (22¾ x 6¾ x 16"H)	49.95

For Imhoff Cones See Page 54. The Complete Line of Scienceware® Funnels Begins on Page 51



Separatory Funnel Rack

Chemical Resistant and Fully Adjustable

This high-density polyethylene rack holds four separatory funnels up to 98mm (3¾") in diameter on its center shelf. Four beakers fit on the bottom shelf and four long-stem funnels on the top shelf. Set screws in the base risers allow vertical adjustment of telescoping rods, while set screws in the center shelf permit adjustment for accommodating separatory funnels of various lengths. Fully extended, the rack measures 203 x 489 x 508mm H

(8 x 19¼ x 20"H). Four feet extend from the base for stability and have holes for permanent mounting.

HDPE

CATALOG NO.	PRICE/EACH
H18967-0000	\$360.00




Rotary Separatory Funnel Racks

Fully Adjustable Support Disc

Molded of strong, rigid white melamine, the rack consists of a circular base, screw-in column and a support disc held in place by a screw clamp. The support disc is fully adjustable to any desired height. Support holes for funnels are open front type for easy removal and placement of separatory funnels. Shipped unassembled.

CATALOG NO.	FOR FUNNELS	PLACES	HEIGHT X DIAMETER	PRICE/EACH
H18963-0000	100cc	10	530 x 300 mm (21 x 11¾")	\$154.50
H18964-0000	500/1000cc	6	530 x 350 mm (21 x 13¾")	135.00

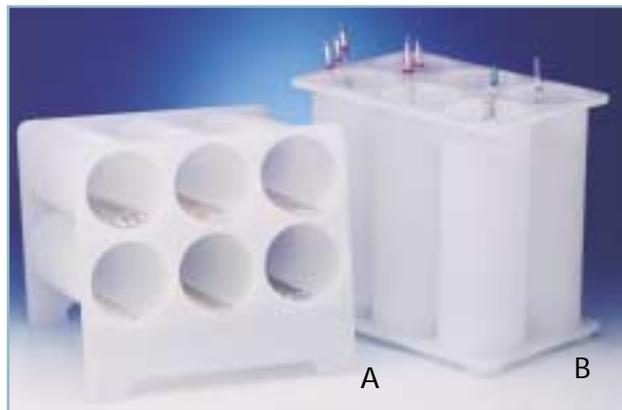

Volume Pipette Storage Racks

Two Models for Vertical or Horizontal Pipette Storage

Heavy-duty polyethylene racks have six tubes that hold 40 to 50 pipettes each. A polyethylene foam disk in the bottom of each 102mm (4") I.D. tube cushions the pipettes and protects tips from damage. The 305mm (12") tube depth stores pipettes safely and provides easy withdrawal. Especially convenient for holding pipettes in a sterile hood environment.

PE

A. The Horizontal Rack is tilted slightly back so pipettes will not slide forward. This rack is 343 x 330 x 318 mm H (13½ x 13 x 12½"H).

CATALOG NO.	PRICE/EACH
H18950-0000	\$387.00

B. The Vertical Rack stands on four feet. The base plate measures 241 x 356 mm (9½ x 14") and the total height is 33cm (13").

CATALOG NO.	PRICE/EACH
F18950-0001	\$360.00

Hydrometer Racks

Safe Storage for Glass Hydrometers

Protect fragile glass hydrometers by storing them in this 18-place polypropylene rack. The middle and top plates have 25mm (1") diameter holes. The bottom plate is 44mm (1¾") high and has wells tapering from 25mm to 9mm (1 to ¾") to hold the hydrometer bulb. The base is 210 x 216mm (8¼ x 8½") for tip-proof stability. Overall height 251mm (9⅞") with handles. Steam autoclavable at 121°C (250°F). *Per each, 6 per case.*

PP A



CATALOG NO.	DESCRIPTION	PRICE/EACH
F18988-0000	For Long Hydrometers Top plate is 178mm H (7"H)	\$57.50
F18988-0001	For Short Hydrometers Top plate is 143mm H (5½"H)	57.50

Thermometer Rack

Safely Holds Up to 25 Thermometers

This polypropylene rack holds up to 25 thermometers in a safe vertical position to avoid damage while drying. Hole diameters in the top and middle plates are 7.9mm (⅝") while the bottom plate has only 3.2mm (⅛") diameter drain holes. The distance between top and bottom plates is 105mm (4⅛"). Overall dimensions are 149 x 213 x 251 mm H (5⅝ x 8⅝ x 9⅞"H). Steam autoclavable at 121°C (250°F). *Per each, 6 per case.*

PP A



CATALOG NO.	PRICE/EACH
F18981-0001	\$32.90